



230 Universal Procedures Table



The Right Thing to Do - For Your Practice, Your Staff and Your Patients

As the practice of medicine changes, more procedures are being performed in the office than ever before, making efficiency essential.

Positioning patients can be a difficult task when dealing with less mobile, elderly and special needs patients. The Ritter 230, with its low height and multipositioning capabilities is designed to improve office efficiency by allowing you and your staff to easily position patients with the touch of a button; reducing the likelihood of fatigue, back and neck strain, and workplace injury.

The Ritter 230 is designed to make patients feel more comfortable and secure when they can easily get onto a table and experience smooth and steady movement during positioning. This allows patients' anxiety to be reduced, and helps them to feel satisfied with the level of care you provide.

Consider the Ritter 230 Power Procedures Table a long-term investment for your practice, your staff and your patients.

The Ritter 230 with power height, back, foot and tilt allows you to effortlessly and precisely position patients for a multitude of procedures.



*Product shown with electrical outlet option & Dusty Blue upholstery

Ergonomic and intuitive standard hand control in easy-to-find location.

*Shown with Moss upholstery

230 Ritter



The Ritter 230 – designed to satisfy your positioning needs.



Knee-Chest Tool-less reconfigure of the foot support provides the perfect position for lower body exams and procedures.



Lithotomy Extending the stirrups and treatment pan, lowering the foot support and tilting the seat provide the perfect position for lower body exams and procedures.



Flat The table supports exams and procedures requiring supine, prone, left lateral or right lateral positioning.

A power procedures table is a wise investment for your practice.

facing you and your practice today, so choosing the right equipment is more important than ever.

More and more healthcare providers are investing in power tables when updating or constructing new facilities. This new trend is making exam and procedure rooms more efficient and safer, giving the patient a better overall healthcare experience.

Make the smart choice. Choose the Ritter 230 Power Procedures Table and experience the many benefits of a power table.



Low Entry With the increasing number of elderly and disabled patients, the low entry position is designed to make it easy and safe for the patient to get on and off the table.



Trendelenburg Whether better positioning for the patient, treating for shock or faintness, or increasing anesthesia onset, the Trendelenburg position offers a full 30° of tilt.



The Ritter 230 is a power procedures table with basic features, making it very affordable. Better yet, the Ritter brand has provided quality to healthcare professionals for over 50 years and is recognized as an industry leader.

*Shown with Clay upholstery

Options designed to improve efficiency





Rotation

Exam and procedure rooms are being equipped with more and larger devices. Inefficiency results when moving these devices to the patient. Outfitting your procedures table with the rotation option simplifies this task by allowing you to move the patient instead of the device.

The simple touch of a pedal (highlighted in photo) on either side of table permits 350 degrees of rotation; touch the pedal again to lock table in place. It's that easy!



A 24-inch wide and flat upholstery option offers greater access to the patient. The seamless design ensures a nice look and easy-to-clean surface.

Flat Upholstery

*Shown with Navy upholstery



Electrical Outlets Optional hospital grade outlets (on each side) allow easy access to electrical power, providing a clean look to the room and floor space free of power cords.

Accessories that add versatility

of tools.

Adjustable Arm System

Designed to offer patients a stable

platform in accessing the table and

sense of security during positioning.

Arms rotate and lock into one of

three possible positions and are

removable, all without the use



Articulating Armboard Designed to comfortably support the patient's arm. Ball joint design permits multi-directional movement.



Armboard Supports patient's arm. The delicate procedures.




24" Wide Foot Rest Designed to provide extra width for added support under the patient's feet.

Articulating Knee Crutch

for lengthier procedures.



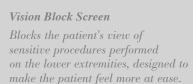




Provides extra side-to-side security for the patient when side panels are in the vertical position. When in the flat position, the side panels can serve as an extra work surface.




Provides mobility without sacrificing valuable floor space.

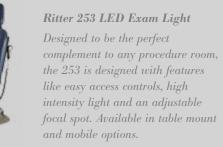






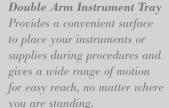



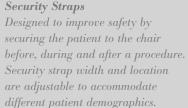

Security Side Panels





telescoping pole rests on the floor to provide extra stability for

Ritter 230

S P E C I F I C A T I O N S



Ritter 230 Universal Power Procedures Table (shown with optional adjustable arm system)

Specifications

- Patient Weight Capacity: 450 lb (204 kg)
- Height: Minimum 22.5" (57.2 cm) Maximum 40" (101.6 cm)
- Seat Tilt Angle Range: 0° to 30°
- Back Support Angle Range: 0° to 85°

Foot Support Angle Range: 0° to 90°

- Upholstery Width Dimensions: Headrest: 18" (45.7 cm) Back support: 28" (71.1 cm) Seat: 28" (71.1 cm) Foot support: 18" (45.7 cm)
- Patient Support Surface Length: Flat position: 72" (185.4 cm) Flat position with headrest extended: 83" (213.4 cm)

Paper Roll Holder Capacity: Holds 18" x 3.5" paper roll (45.7 x 8.9 cm)

Stainless Steel Treatment Pan Dimensions: 11.5" L x 9" W x 2.5" D (29.2 x 22.9 x 6.4 cm)

Stirrups are standard

Optional Electrical Outlet Rating: 115 VAC, 3 Amps

Accessories

Ritter 253 LED Exam Light Foot control - linear style

24" wide foot support

12" wide headrest

Adjustable arm system (right and left side arms included)

Adjustable arm receptacle

Articulating armboard

Hand surgery armboard Fixed armboard

Knee crutches - articulating and fixed Casters¹

Caster base^{1, 2}

Security straps

Security side panels Base rails

Instrument tray

IV pole Urology drain system¹ Special procedures headrest Facial pad Vision block screen

3' power cord

Options

Electrical Outlets: Duplex, hospital grade¹ Rotation: Models equipped with rotation option add 3 1/2 inches to minimum and maximum height specifications

Flat upholstery set, 24" wide

Electrical Ratings

115 VAC, 10 Amps, 50/60 Hz (for electrical outlet models)

115/230 VAC, 7/3.5 Amps, 50/60 Hz (for non-electrical outlet models)

115/230 VAC, 7.5/3.8 Amps, 50/60 Hz (for rotation models)

Safety Standards

UL 60601-1 and CAN/CSA C22.2 No. 601.1 - M90

Ritter Upholstery Colors:

Navy-231 Dusty Blue-233 Moss-230 Clay-234 Pebble Grey-216 Shadow-232

Perfect Plum-229

Black-312

¹Not available on rotation models ²Not compatible with 24" wide foot support

CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company.

For more information or a demonstration, contact your Ritter or Midmark Dealer, or contact Midmark Corporation, Versailles, Ohio. Call: 1-800-MIDMARK Fax: 1-800-365-8631 midmark.com

© Midmark Corporation 2004 Products subject to improvement changes without notice. Litho in U.S.A. 007-0544-00 Rev. H1 (9/14)

